



Traceable® Thermohygrometer with Clock

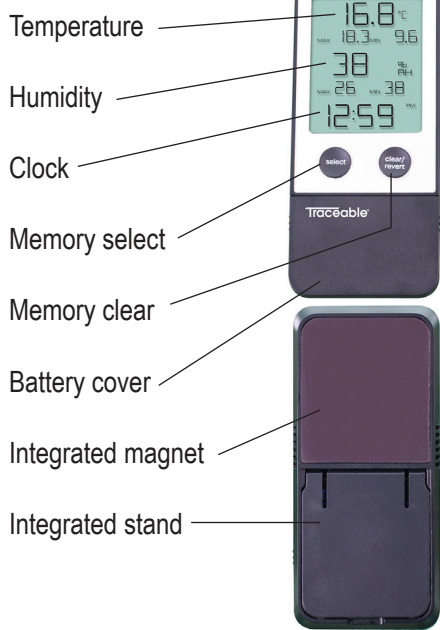


INITIAL SETUP


- 1** Remove battery cover and insert batteries.

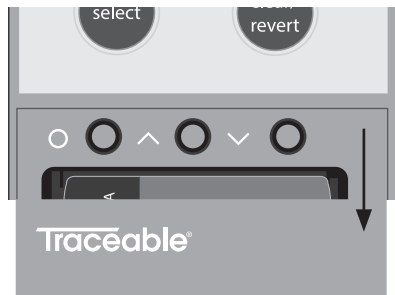

- 2** Once batteries are installed, press **O** button to select setting.

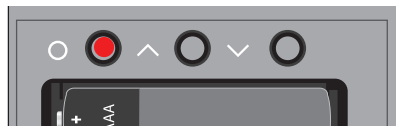


- 3** The first setting, Hours, will flash on display.

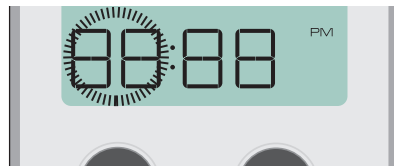

- 4** Press “^” and “v” buttons to adjust setting.




- 5** Once finished, press **O** button to save setting and proceed to Minutes, 12/24 Hour, etc.

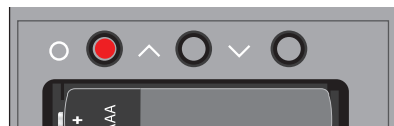


- 6** Complete all settings to configure device:

- Hours
- Minutes
- 12 or 24 hour format clock
- Temperature units

- 7** Replace battery cover.

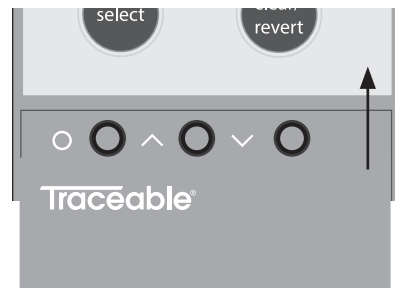

DISPLAY DETAILS

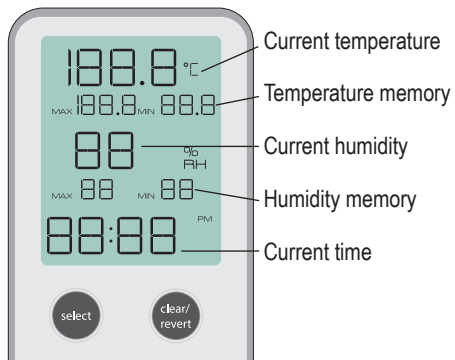

CLEARING INDIVIDUAL READINGS

1. Press select button until desired reading to be cleared is flashing.
2. Press and hold clear/revert button for 3 seconds to clear reading.

CLEARING ALL MEMORY READINGS

1. Press and hold clear/revert button for 3 seconds to clear all memory readings.

WARRANTY, SERVICE, OR RECALIBRATION

Traceable[®]

an antylia scientific company

12554 Old Galveston Rd. Suite B230
Webster, Texas 77598 USA
Ph. 281-482-1714 • Fax 281-482-9448
support@traceable.com
traceable.com

Traceable[®] Products are ISO 9001:2018 Quality-Certified by DNV and ISO/IEC 17025:2017 accredited as a Calibration Laboratory by A2LA.

Item No. 56000-10

Model No. 5660

©2022 1065T6_M_92-5660-00 Rev. 0 09022022

Traceable[®] Thermohygrometer with Clock Quick-Start Guide

